

S-CAPEPOD Evacuation Sheet

The S-CAPEPOD Evacuation Sheet is a swift, life saving evacuation aid for non-ambulatory patients that is permanently placed underneath a patient's mattress. The patient neither sees it, nor feels it. It is intuitively designed so the patient will be ready for evacuation in a matter of seconds with a few simple steps. Select standard or bariatric.



Features

- Allows 1-4 people (depending on standard or bariatric model) to evacuate non-ambulatory people along with the bed mattress. Use this life-saving tool to minimize the likelihood of leaving people behind when disaster strikes
- A unique packing method helps rescuers to act quickly. The S-CAPEPOD rescuers wrap the patient on their mattress by means of three body flaps. As a result, the patient lies in a protected cocoon while being evacuated
- S-CAPEPODs are fire retardant and anti-bacterial
- S-CAPEPOD permanently remains underneath bed mattress and ready to use at all times
- Bariatric model has an extra abdominal belt that firmly fixates the patient onto the mattress to prevent patient from moving during a vertical evacuation such as down a staircase
- Suitable for all different types of mattresses
- Easy to clean with a disinfectant wipe, or you can place the S-CAPEPOD in a washing machine to clean
- User instructions are printed on flaps
- Ideal for non-ambulatory patients and people with reduced mobility

Evacuate Non-ambulatory Patients Quickly and Safely in a Few Easy Steps



Simply release the flaps of the S-CAPEPOD from beneath the bed, and fold over patient as directed in written manual. Then slide the patient and mattress off the bed using the S-CAPEPOD towing handles, and evacuate to safety.


See next page for specs >

S-CAPEPOD Evacuation Sheet

Where to place the S-CAPEPOD

Place the S-CAPEPOD directly on top of the bed frame so it is always ready to use. The S-CAPEPOD lays permanently between the bedframe and mattress. The patient doesn't see it nor feels it.

Then, attach the S-CAPEPOD to the bed mattress using the strap adjusters at the top end of the S-CAPEPOD.

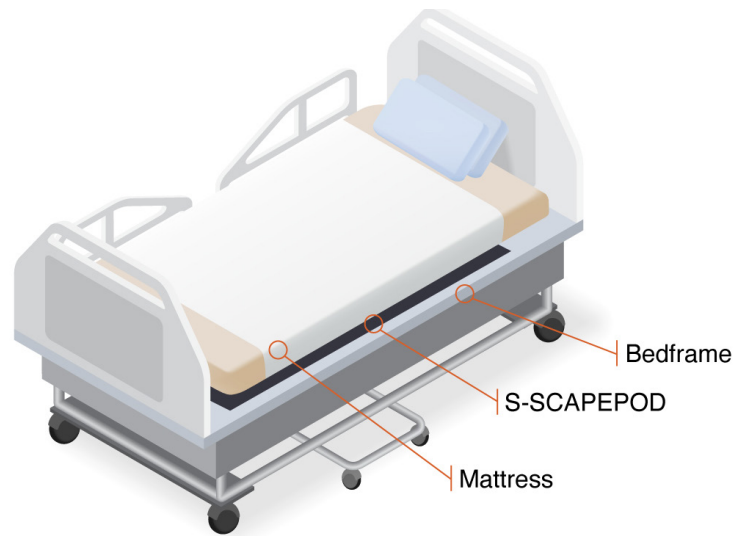

	S-CAPEPOD Standard	S-CAPEPOD Bariatric
Size under mattress	82" x 33"	82" x 33"
Suitable for mattress sizes	Width: 31" to 38" Maximum length: 82"	Width: 39" to 49" Maximum length: 85"
Safe working load	264 - 330 lbs depending on body size	660 lbs
Product weight	2.425 lbs	5.5 lbs
Material	100% polyester (fire retardant & anti - bacterial)	100% polyester (fire retardant & anti - bacterial)
Classifications	Class 1 according to Regulation MDR 2017/745/EU. Applied standards: NEN-EN-ISO 13485:2016 EN NEN-EN-ISO 14971:2019 EN	Class 1 according to Regulation MDR 2017/745/EU. Applied standards: NEN-EN-ISO 13485:2016 EN NEN-EN-ISO 14971:2019 EN
Users needed	1-2 people	1-4 people
Guarantee	1-year warranty against defects in materials and production, which is valid from the date of purchase. Improper use, improper maintenance, vandalism, product modifications or accidents are not covered by the warranty.	1-year warranty against defects in materials and production, which is valid from the date of purchase. Improper use, improper maintenance, vandalism, product modifications or accidents are not covered by the warranty.

Cleaning and Maintenance

By hand: Use a wipe and a disinfectant

In the washing machine: Washing temperature, 40 degrees (Celsius), dry at a low temperature, do not iron, do not dry clean, do not bleach, do not use washing powder with fabric softener.